



THANK YOU

Before I tell you how Easy Smile works I really want to thank you for the purchase you have just made.

We at the Upshell company work hard every day to solve your needs We try to offer the best products at a price everyone can afford. If you find Easy Smile useful and effective after the first application, we ask you a small favour: Share your result via a review on Amazon: you will help other people us to grow as an Italian company, but also other people who are looking for the best whitener among those on sale on Amazon.



HOW IT WORKS?

Our Easy Smile whitener adopts a special innovative blue LED light technology. The blue rays will react with the whitening gel and effectively remove stains on the surface of the teeth. The mouthpiece has been designed through engineering studies to be as soft as possible. You will find the Easy Smile Kit comfortable and feel no discomfort during the whitening process. When you bite down on the mouthpiece, it fits and completely covers your teeth. Before starting the treatment, it is very important to know how it works and what steps to follow. Please read the instructions carefully before starting the treatment!

APPLICATIONS

The bleacher is ideal for the removal of: Coffee stains, Tobacco, genetic stains, tetracycline, food dyes, ageing and other factors that can cause discolouration.

STORAGE

Do not expose it to direct sunlight and store it in a cool place.

USER GUIDE

1. Rinse the mouthpiece and mouth before use.
2. Look in the mirror to compare your teeth with the guide paper and check the starting colour. 3. Rotate the teeth whitening gel pen until the gel flows towards the brush. Apply the gel evenly to all teeth, both upper and lower. Do not use too much gel, more gel will not give a better or faster result. Use about 0.5 ml - 0.7 ml each time. Avoid applying the gel directly on the gums or inside the lips. If you mistakenly apply it on the gums, do not be alarmed as the burning sensation will disappear within a few minutes.

Be sure to rinse with water.

4. Make sure the charge is complete before first use. Put the mouthpiece in your mouth and bite down on it, then switch on the LED light by holding down the button. Breathe deeply to reduce saliva production. (Excessive saliva dilutes the gel concentration).

5. The whitening device switches off automatically after 16 minutes. Take it out, rinse the mouthpiece and mouth with water and compare your teeth again with the colour guide card to see the improvement.


6. Do not eat or drink anything, except water, for 30 minutes after finishing the treatment.



Rinse the mouth
before use



Apply the gel with the pen
on all teeth

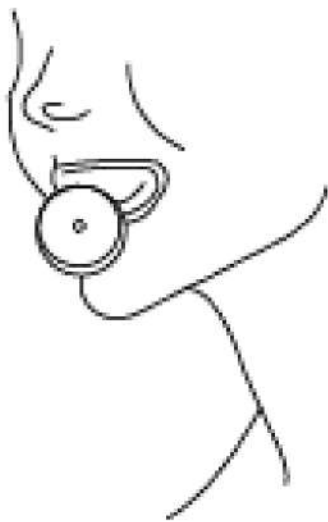

Insert the mouthpiece and bite
down on it, then switch on the LED
light



Wait 16 minutes to
conclude the treatment



Rinse the mouth
with water



Wash mouthpiece with water
and close gel pen with cap

REMINDER

1. Rinse the mouth thoroughly, completely dry the water on the surface of the teeth and apply the gel to the teeth evenly.
2. Excess saliva may dilute the concentration of the gel.
3. Everyone's teeth are different, some people experience slight sensitivity, it is a normal phenomenon. A desensitising gel can be applied to relieve discomfort. If the discomfort is severe, reduce the time of use or discontinue use immediately.
4. Avoid applying the whitening gel directly to the gum area or inside the lips.
5. If the gel inadvertently gets on your face or clothes, wash it off immediately.
6. Please put the cap on the top of the pen after finishing the teeth whitening treatment to prevent the active ingredient from volatilising.
7. Store the product in a place out of the reach of children.
8. It is recommended not to use the whitener more than 4 times a week

INGREDIENTS

Sorbitol, Water, Phthalimidoperoxyproic acid (PAP), Propylene glycol, Glycerine, Potassium nitrate, PEG-8, Hydroxyapatite, Cellulose gum, Hydroxyethylcellulose, Xanthan gum, Sodium saccharin, Menthol, Methylparaben, Aloe Barbadensis (Aloe Vera) leaf juice, Chamomilla Recutita extract, Punica Granatum seed extract, Sodium bicarbonate, Propylparaben.

NOTES

1. Charging time is about 2 hours.
2. Charging status: The Easy smile mouthpiece will flash when in charging status
3. Press the button 2 seconds to switch on/off.
4. The service life of the whitener is approximately 24 months